

ice point™ Calibration Reference Chamber

High Precision Thermoelectric Chamber Uses Water to Maintain 0°C (32°F) Continuously

TRCIII-A



Standard 230 Vac model only

- ✓ Stylish Design
- ✓ Portable Rugged Case for Use in the Factory, Calibration Laboratory, or Instrument Shop
- ✓ Built-in Temperature Readout with 0.1 (°F/°C) Resolution
- ✓ Accepts up to 6 Probes at One Time
- ✓ Calibration of All Probe Type Thermometers and Sensors



TRCIII-A shown smaller than actual size.

Probes sold separately;
See next page.

The ice point™ TRCIII-A Calibration Reference Chamber is part of OMEGA's fine line of calibration reference instrumentation.

The TRCIII-A ice point™ Reference Chamber relies on the equilibrium of ice and distilled, deionized water at atmospheric pressure to maintain 6 reference wells at precisely 0°C (32°F). Each well extends into a sealed cylindrical chamber containing the distilled deionized water. The outer walls of the chamber are cooled by thermoelectric cooling elements. The increase in volume produced by the creation of ice crystals within the cell is sensed by the expansion of a bellows, which operates a microswitch and controls the cooling elements. The alternate freezing and thawing of the ice accurately maintains a 0°C (32°F) environment around the reference wells.

Features

- ✓ 230 Vac CE Marked Models
- ✓ NIST Calibration Certificate Included
- ✓ Bench-Top Design with Tilt Handle Standard
- ✓ Illuminated Power Switch
- ✓ Companion hot point® Dry Block Probe Calibrator Available (Model CL900A/950A)

Specifications

Reference Wells: 6 wells, 4.0 mm ($\frac{1}{8}$ " ID, 95.25 mm ($3\frac{3}{4}$ " deep; wells are thermally and electrically grounded to each other; accepts 3 mm ($\frac{1}{8}$ " OD probes

Reference Temperature: 0°C (32°F)

Accuracy: $\pm 0.1^\circ\text{C}$ (± 0.18)

Stability: $\pm 0.03^\circ\text{C}$ (± 0.07) for constant ambient

Environmental Ambient:

2 to 32°C (35 to 90°F)

Power: 115/230 Vac, 50 Hz, 104 W

Stabilization Time: 2 hours after initial power up

Dimensions:

203 H x 419 W x 305 mm D
(8 x 16½ x 12")


Weight: 11.8 kg (26 lb)



OMEGACARESM extended warranty program is available for models shown on this page. Ask your sales representative for full details when placing an order. OMEGACARESM covers parts, labor and equivalent loaners.

Caution:
Unit must not be allowed to freeze or damage will result.

We make running changes when technical advances allow. Check at time of ordering for additional features.

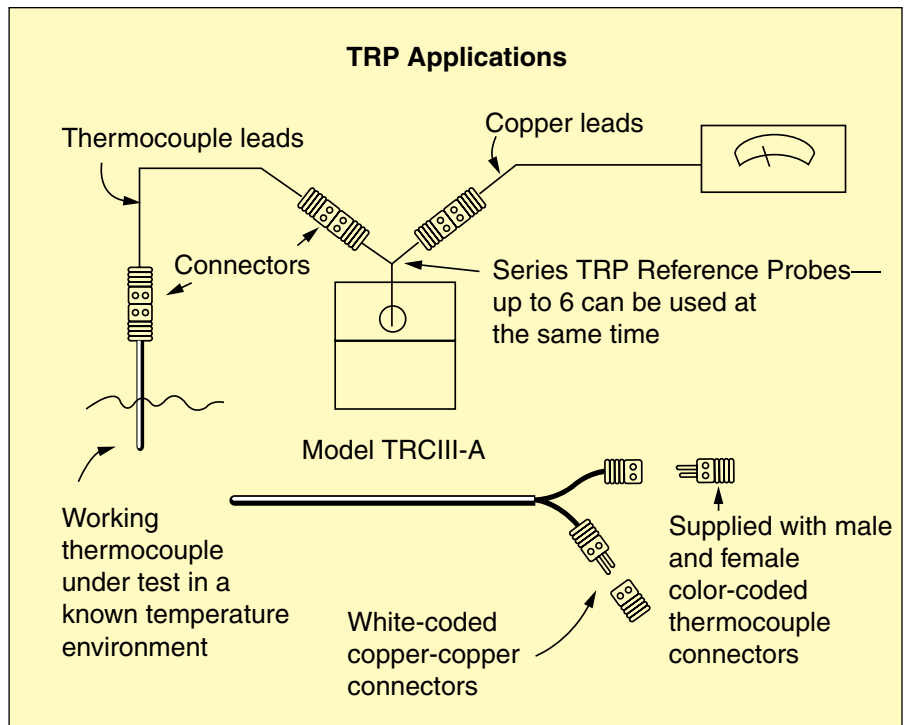

To Order

Model No.	Description
TRCIII-A	ice point™ reference, 120 Vac operation
TRCIII-A-230VAC*	ice point™ reference, 230 Vac operation


Comes complete with operator's manual and NIST calibration certificate.

* **Note:** Only 230 Vac model is available as CE marked.

Ordering Examples: TRCIII-A, ice point™ reference cell, TRP-K, TRP-J, TRP-T, TRP-E, Type K, J, T and E temperature reference probes.

OCW-3, OMEGACARESM extends standard 1-year warranty to a total of 4 years.

TRCIII-A Accessories

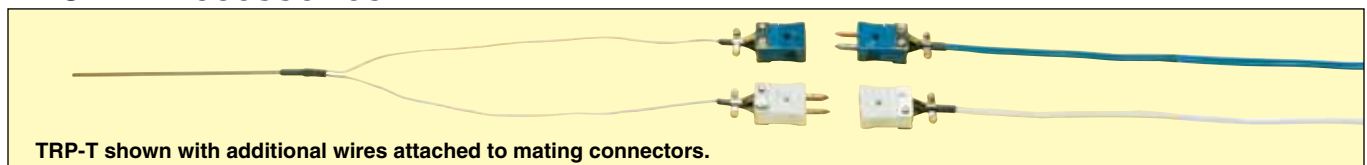

Thermocouple TRP Features

- ✓ Made with Special Limits of Error Material
- ✓ All Standard Calibrations Available
- ✓ Comes Complete with Mating Connectors
- ✓ Can Be Ordered with 12, 24 or 36" Leads

Thermocouple TRP Probes

Model No.	Thermocouple Material Type	
	+	-
TRP-K	CHROMEGA®	ALOMEGA®
TRP-T	Copper	Constantan
TRP-E	CHROMEGA®	Constantan
TRP-J	Iron	Constantan
TRP-S	Platinum	Platinum 10% Rh
TRP-R	Platinum	Platinum 13% Rh
TRP-B	Platinum 30% Rh	Platinum 6% Rh
TRP-C	W 5% Re	W 26% Re
TRP-D	W 3% Re	W 25% Re

* Refrigeration required for below ambient cooling.

** Available as an option where not supplied as standard as HCTB-COIL.

†† For 240 Vac power, add suffix "-240V" to model number, no additional charge.

Note: For 4-point NIST calibration, points: -20°C, 60°C, 120°C, 180°C, add suffix "-NIST" to model number for additional cost.